



STEAMBATH

Model: 401



FEATURES

- Laboratory steambaths or hot water-baths are suited for many types of digestion experiments
- The self locating removable lids can readily accommodate a variety of round bottom kjeldahl flasks
- Manufactured from high grade 304 stainless steel
- The unique design of the lid and tank rim ensure that all water vapor returns back to the bath
- Easily replaceable high quality incoloy heating element for long and reliable use
- Manual reset safety thermostat protects the element should the bath run dry
- Simple device for automatic control of water level
- Flask openings 31mm, 50mm, 70mm and 80mm diameter
- Heat resistant, safe to touch lid knobs

NOTE

- 6 station unit



ALL SCIENTIFIC EQUIPMENT COMES WITH A TWO
YEAR WARRANTY ON COMPONENTS AND
DEFECTIVE WORKMANSHIP - Ts & Cs APPLY.

www.scientific.co.za






STEAMBATH

Model: 401

SPECIFICATIONS

POWER	2000W
VOLTAGE	230V / 50Hz
PACKAGING (W x D x H)	550 x 400 x 245 (mm)
UNIT (W x D x H)	545 x 280 x 160 (mm)
NET WEIGHT	5kg
GROSS WEIGHT	6kg
SHIPPING VOLUME	0.05m ³
NUMBER OF STATIONS	6

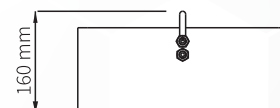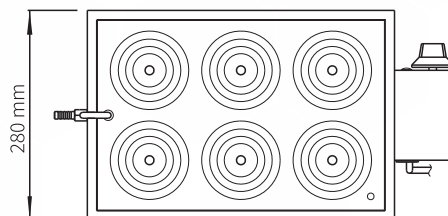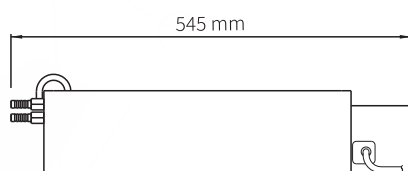

NOTE

- When viewing the appliance from the front in its usual operating position, the width (W) of the product is the total distance from left to right; the depth (D) of the product is the total distance from the front to the back; the height (H) is the total distance from the bottom of the product to the top. Pictures may vary from actual product.

DISCLAIMER

- Output performance figures quoted are dependent on various factors. Scientific reserves the right, without notice, to make changes and revisions to product specifications, materials and design, when we believe it will provide better performance, durability, and efficiency. Unit measurements may vary depending on point of reference.

WARRANTY

- All Scientific equipment comes with a standard ONE year warranty. The 1 Year Extended Warranty is only valid for customers who have registered their warranty online at www.scientific.co.za within the first year after making their purchase.



ALL SCIENTIFIC EQUIPMENT COMES WITH A TWO
YEAR WARRANTY ON COMPONENTS AND
DEFECTIVE WORKMANSHIP - Ts & Cs APPLY.

www.scientific.co.za

Published Date: 2020.09.25